

112C5b (9024) Transducer

The I12C5b intraoperative simultaneous biplane transducer offers detailed imaging of difficult-to-access anatomy, with I and T shaped arrays.

- Biplane transducer with simultaneous imaging in two planes for precise guidance
- Compact size for easier access to tight, difficult-to-reach sections of the abdominal cavity. Switch between arrays to capture the clearest image
- Flexibility to take freehand biopsies
- Advanced transducer connector designed to reduce noise in the ultrasound image, and for convenient one-handed operation
- Perform quick, reliable puncture and biopsy procedures with our single-use needle guide
- Unique, fully-immersible transducer for easy sterilization



Fig. 1 I12C5b transducer

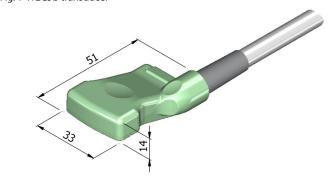


Fig. 2 Technical drawing and dimensions (mm)

Specifications

Transducer Name		I12C5b
Transducer Reference Type		9024
Compatible Ultrasound Systems		bk3500, bk5000
Applications (Typical)		Intraoperative ¹ , Musculoskeletal (conventional), Peripheral Vessel (peripheral vascular)
Basic Imaging Modes ²		B, M, CFM, Doppler, Tissue Harmonic Imaging, Contrast Imaging, Elastography
Imaging Planes		Transverse and Sagittal
Weight (Approximate)	g	45
Acoustic Aperture	mm	23 x 5
Number of Elements		2 x 96
Radius of Curvature	mm	40
Image Plane Focal Length	mm	Variable
Image Field		Sector 32° Expanded Sector 55.5°
Focal Range	mm	10 - 80
Frame Rate (Max)	Hz	>200

¹ bk5000 only. Liver, biliary, pancreas, kidney. 2 Not all systems support all imaging modes.

Specifications

SAFFTY

Designed and tested in accordance with current version of EN60601-1 (IEC60601-1) "Medical Electrical Equipment, Part 1". When used with BK Medical's systems, Type BF requirements are met.

ENVIRONMENTAL

Operating Pressure:

700-1060 hPa (normal atmospheric pressure)

Operating Temperature:

+10 to +40 °C (+50 to +104 °F)

Storage Temperature:

-25 to +70 °C (-13 to +158 °F)

REPROCESSING

- The reprocessing temperature should not exceed +60 °C.
- Complete details and procedures can be found in *Care and Cleaning*.
- Follow manufacturer's instructions for use.

VALIDATED CLEANING AND DISINFECTION

For all cleaning and disinfection methods except wiping, plug must be covered.

Manual Cleaning

3E-Zyme™

Manual Disinfection

Ethanol 70% (wiping) Korsolex® Basic

Revital-Ox™ RESERT®/RESERT® XL HLD

Tristel® Fuse for Stella

Automatic Disinfection¹

MEDIVATORS® Advantage® Plus, with:

- MEDIVATORS Intercept® Detergent
 MEDIVATORS Rapicide® PA
- 70% Isopropanol Flush

VALIDATED STERILIZATION¹

Immersion Methods (plug must be covered)
STERIS® SYSTEM 1®, SYSTEM 1E®, SYSTEM 1 Plus,
SYSTEM 1 Express

Gas Processes (plug must be uncovered)

STERIS® V-Pro™
STERRAD NX™ and 100NX Systems
STERRAD 100S and 200 Systems

POWER SUPPLY

Internally from system

CABLE LENGTH

2.8 m (9.2 ft)

¹ Automated reprocessing methods validated for 100

Ordering Information

ACCESSORIES AVAILABLE

KE4300: Carrying Case

UA1414: Leakage Testing Kit

UA1337: Single-use needle

Single-use needle guide. Sterile packed. Includes 2 needle guide insert palettes (3-angle and free-angle respectively), each with 9 different inserts for 10-20

gauge needles (pack of 9)

UA2435: Reprocessing Tray

TRANSDUCER COVERS

UA0004¹: Intraoperative sterile cable covers, latex-free, 7.6 x 244 cm (pack of 24)

UA0008¹: CIV-Flex™ latex-free covers.

sterile, 11.9 cm tapered to 3.8 x 61cm, (pack of 24). Includes gel packet

UA0073: Sterile, latex-free, 15.2 x 244 cm,

telescopically-folded with attached 3

cm NeoGuard™ tip (pack of 12)

¹Must be used together.


Accessories

Single-Use Needle Guide UA1337



- **Dimensions:** 35 x 22 x 10 mm (1.4 x 0.9 x 0.4 in)
- Weight: 4 g
- Materials: MABS plastic with 20% Barium sulphate

Single-Use Needle Guide Kit UA1337 – Needle Guide Insert Palettes

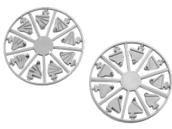

- **Dimensions:** Ø120 X 11 mm (Ø4.7 x 0.4 in)
- Weight: 28 g
- Material: MABS plastic with 20% Barium sulphate

BK Medical ApS, Mileparken 34, 2730 Herlev, Denmark.T +45 4452 8100 F +45 4452 8199

BK Medical 8 Centennial Drive Peabody MA01960 T + 1 978-326-1300 bkmedical.com

USA Sales & Service BK Medical 8 Centennial Drive Peabody MA01960 USA T + 1 978-326-1300 F + 1 978-326-1399

bkmedical.com

Europe and Rest of World Sales, Service & Design Center BK Medical Mileparken 34 2730 Herlev Denmark T +45 4452 8100 F +45 4452 8199 bkmedical.com

China Sales & Service Analogic Medical Equipment (Shanghai) Co., Ltd. 1377, Lan Dian Road Pu Dong New District Shanghai China 201132 T +86 21 2089 0333 bkmedical.com

bk medical •))